

BeneView™ T5

Patient Monitor



Handle
Convenient for portability



Touch Screen
Directly access many functions through simple one-touch commands



CF Storage Card
4 GB storage capability
Stores up to 120 hours of graphic / tabular trend and 48 hours of full disclosure waveform for each patient



Abundant Connectors
For keyboard, mouse, analog output, networking and slave display



Thermal Recorder
Built-in 3 channel recorder offers convenient documentation



Battery
Up to 5.5 hours working time with two rechargeable Li-ion batteries



Wall Mount



Bedrail Clamp



Rolling Stand



Anesthesia Machine Mount



BeneView™ T5

Patient Monitor



DISTRIBUTOR:



Masimo is the trademark owned by Masimo Corporation. Nellcor is the trademark owned by Nellcor Puritan Bennett LLC. MINDRAY™ is a trademark of Shenzhen Mindray Bio-Medical Electronics Co., Ltd. All rights reserved. © 2010 Shenzhen Mindray Bio-Medical Electronics Co., Ltd. All rights reserved. P/N: ENG-BENEVIEW T5-420285-20100125


mindray

Mindray is listed on the NYSE under the symbol "MR"
Mindray Building, Keji 12th Road South, High-tech Industrial Park,
Nanshan, Shenzhen 518057, P.R. China
Tel: +86 755 26582888 Fax: +86 755 26582680
E-mail: intl-market@mindray.com Website: www.mindray.com


mindray
healthcare within reach

BeneView™ T5


Patient Monitor

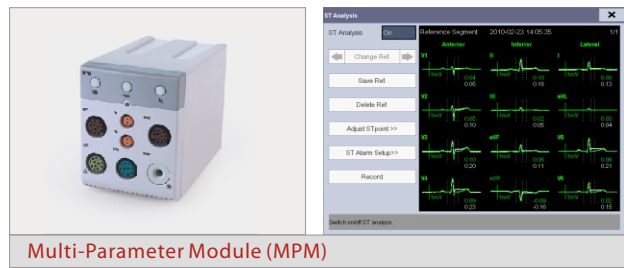


- Lightweight and compact
- Stores up to 24 hours of data
- Conveniently for patient transport and transfer
- Saves space by integrating the most consistently required parameters

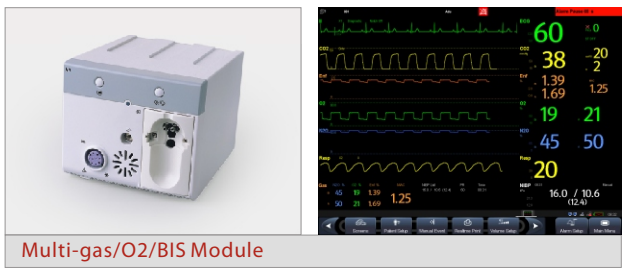


- World-leading technology, smallest size
- Measures and auto identifies all five anaesthetic agents
- MAC (Minimum Alveolar Concentration)
- Agent mixture detection
- Optional integrated BIS module

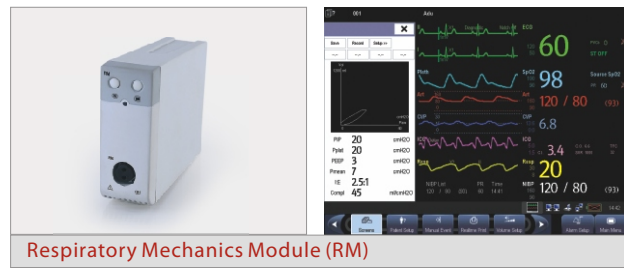

- Provides continuous and noninvasive monitoring for spirometry parameters
- Detects any change in patient's airway
- Assists clinicians to manage ventilation devices

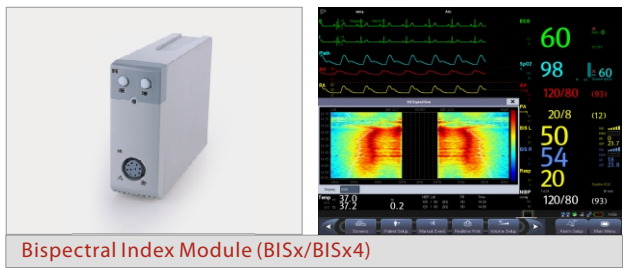


- 4-ch EEG waveforms, DSA
- FDA cleared awareness indication for use
- Faster, more predictable wake-ups, improves patient care
- Reduces the risk of intra-operative awareness with recall
- Reduces costs of anesthetics and sedatives

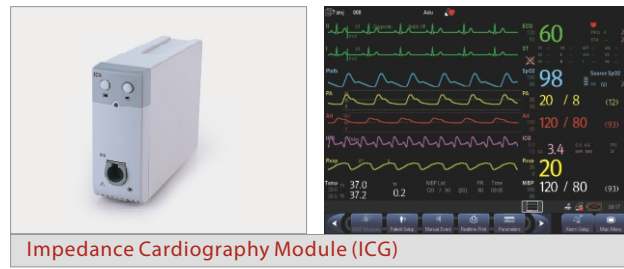


- Offers non-invasive and continuous monitoring of patient's hemodynamics status
- Guides selection and optimization of medications
- Easy to use and very cost effective
- Virtually no patient complications



- Acquires real-time information from Vigilance II or Vigileo
- Provides a wide range of hemodynamic parameters including CCO, SvO₂, SVR, EDV, etc.
- Continuous assessment of CO by assessing DO₂/VO₂ balance with SvO₂ monitoring.
- Continuous assessment of components of stroke volume (preload, afterload, and contractility)

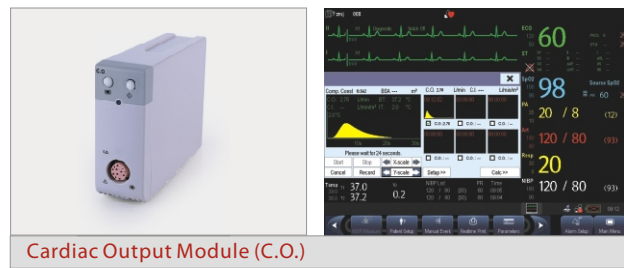

- Traditional right heart thermodilution technology
- Clinical "gold standard" for determining cardiac output
- Displays thermodilution curve and cardiac output, cardiac index, and blood temperature as numerics
- Up to six C.O. measurements can be edited for an average C.O.

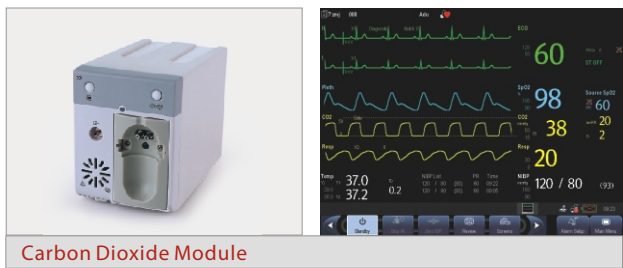

- Mainstream, Microstream or Sidestream EtCO₂ technology for different clinical demands
- Obtains fast, simple and accurate CO₂ results for intubated or non-intubated patients

- 12.1" color TFT display with optional touchscreen
- 13 module slots for flexible configuration
- Advanced technologies including RM, ICG, BISx/BISx4h, CCO, SvO₂, etc.
- Maximum 8 channels for IBP measurement
- Robust data storage capabilities including:
 - 120-hour graphic and tabular trends of all parameters
 - 48-hour full disclosure waveforms review
 - 120 seconds frozen waveforms
 - 1000 pieces of NIBP record storage
 - 100 pieces of regular alarms review
 - 100 pieces of ARR alarms storage
- Hard-wired or wireless network connectivity
- Private mode, night mode, CPB mode
- Transfer configuration through USB dongle



Trend Review
Maximum 120-hour graphic and tabular trends



Big Fond Display
Selectable parameters for distance view



12 Lead ECG
Simultaneous real-time display of 12 leads and up to 12 lead ST segment