


QA3™

Ophthalmic Stretcher

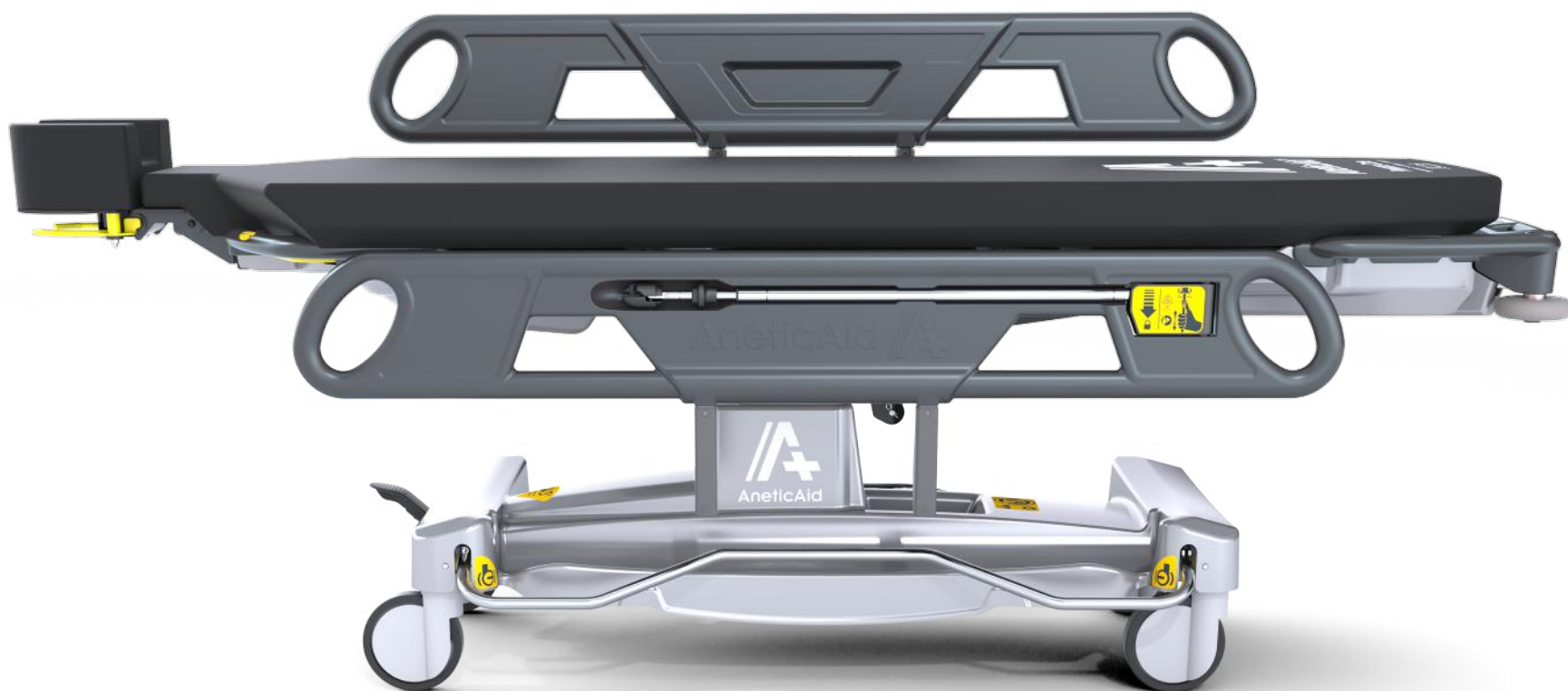

Outstanding surgical access –
exceptional patient comfort



AneticAid® 

STANDARD FEATURES

Spacious leg room, optimal positioning, optimal access from superior or temporal approach – the QA3 Ophthalmic Stretcher offers both the surgeon and the patient an improved surgical journey.

Dual-articulating headrest

This innovative, intuitive and easy to use head positioner is now fitted as standard, articulating not only on the head but the neck as well. The positioner offers superb surgical positioning, surgical access and patient comfort.

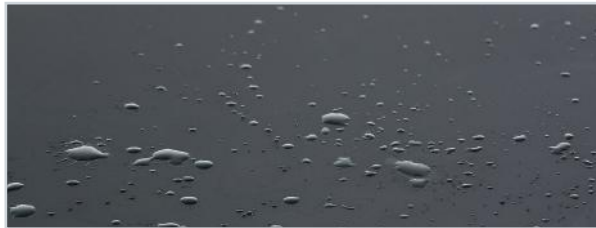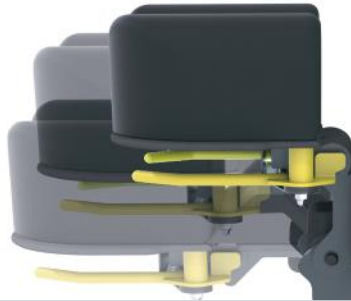

Infection control

Smooth moulded surfaces and acrylic capped vac formings enclose trolley mechanics to minimise dirt traps. Anti-bacterial properties and innovative design for easy cleaning.

And for support workers it contributes to the drive in reducing manual handling with its exceptional low height, powered functions and lightweight manoeuvrability.

Exceptional surgical access

Providing 585mm of clearance underneath the head & back section for comfortable leg room and surgical access from both a superior and temporal approach.

With a low-profile depth sculpted backrest, sculpted shoulder area free from function operating handles and it's low-profile positioning headrest, exceptional surgical access and leg room are achieved.

Press on / Press off 5th wheel steer

Our unique mechanism ensures maximum manoeuvrability and means the QA3 is easily guided by one support worker.

Simultaneous braking low rolling resistance castors

Polyurethane tyres offer low rolling resistance for high manoeuvrability: 150mm diameter central total lock castors fitted as standard, helping to offer a rigid and reliable work platform.

Vertically operated safety side rails

Store away under the mattress improving the ease of mounting and dismounting the stretcher for the patient.

Exceptional patient comfort

Low height of 495mm (excluding mattress) & store away side rails allow for easier mount/dismount, particularly for frail and elderly patients. Also including a versatile positioning headrest, wide patient platform and pressure care mattress for comfort.

Moving and handling

Zero transfers mean risk is minimised for patient and practitioner alike: from the ward, through anaesthetics, to the surgical procedures, and finally to first and second recovery.

Fixed transfusion pole


Innovative quick release safety locking mechanism, integrated into the right-hand foot end side rail for easy access and storage.

Cylinder trough

Side loading for ease of access, suitable for wide range of cylinder options including E, F, CD and HX.

Space saving design

With the back rest and KneeFlex raised, the stretcher has an incredibly small footprint - making it a comfortable and manoeuvrable recovery or transport position and easier to store when not in use.

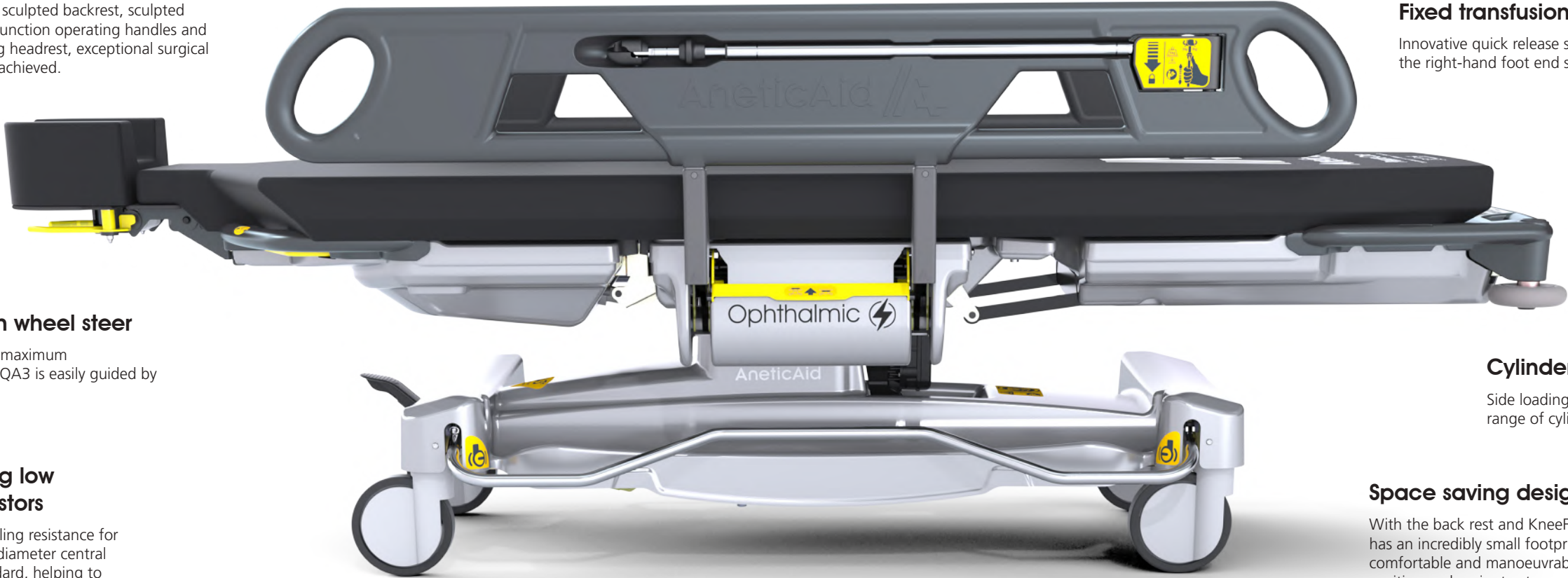

Pressure care mattress

Manufactured from temperature-sensitive gel-feel foam that moulds to the patient's body shape, relieving pressure by distributing the load as evenly as possible. Outer material is anti-microbial, water resistant, vapour permeable and latex free. Now 5% resistant to chlorine and other disinfectant compounds.



AS STANDARD

SPECIFICATIONS		Code: 21118	
Minimum height	*495mm	Mattress length	2245mm
Maximum height	*875mm	Mattress width	700mm
*excludes mattress		Mattress depth	90mm
Surgical Access	585mm	Maximum patient weight	320kg
Minimum length	1700mm	Backrest	0 to +70°
Maximum length	2275mm	KneeFlex	0 to +30°
Overall width	805mm	Trendelenburg	0 to -12°
Castor diameter	150mm	Reverse Trendelenburg	0 to -10°
Weight	135.5kg	Headrest	0 to -87 / +45°

A unique combination of manoeuvrability and stability



DESIGN & INNOVATION



MADE IN THE UK



GLOBAL DISTRIBUTION

⚡ POWERED FEATURES

Minimise risks for practitioners dealing with increasing patient weights with powered positioning: the whole stretcher platform or the backrest can be raised or lowered electronically from the panel at the foot end.

With six powered functions – platform raise and lower, back rest, KneeFlex, Profile, Trendelenburg and Reverse Trendelenburg - as well as CPR return – the QA3™ Ophthalmic Stretcher is a truly powerful and versatile platform for transport, treatment and recovery.



Battery powered

Long-life rechargeable lithium-ion battery (NMC) with highly accurate charge indicator



Electrically Powered Functions


Easy to operate using the hand control positive feel push button controls.



Platform height



Offers smooth movement to position the stretcher at any point in its height range – 495mm to 875mm

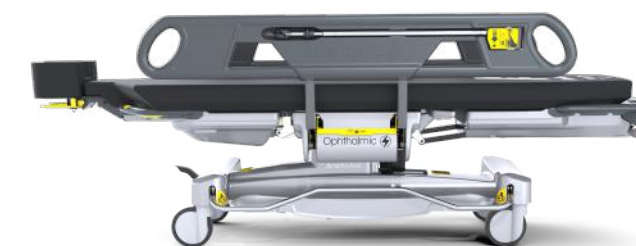

Backrest

Positions the backrest anywhere between the vertical and horizontal for maximum patient comfort



KneeFlex

Easy elevation of the knees to offer improved comfort in a supine position



Profiling

Easy elevation of the back rest and knees at the same time to offer improved comfort during recovery



CPR Return

If a patient needs CPR treatment, there is no time to lose: the powered CPR feature returns the stretcher swiftly to a horizontal, optimal height position for CPR treatment



QA3™ Ophthalmic Stretcher

Foot-end Control

Even greater flexibility when repositioning the QA3™ using the foot-end panel



Trendelenburg / Reverse Trendelenburg

Extremely smooth and controlled by a powered actuator, even at lowest height



DESIGN & INNOVATION



MADE IN THE UK



GLOBAL DISTRIBUTION

Effortless positioning at the touch of a button

ACCESSORISE YOUR QA3 OPHTHALMIC STRETCHER

Customise your QA3 Ophthalmic Stretcher at the time of ordering to better suit your department needs or add greater versatility to your QA3 with this range of convenient accessories.



Oxygen Delivery Bar & Drape Screen

Two products in one. An over patient drape / anaesthetic screen and oxygen delivery system. This versatile flexible accessory mounts onto the side rail via a clamp and can be bent over the patient for drapes to create a screen. The barbed end allows for oxygen tubing to be connected and the fenestrated end delivers a quiet, consistent, flow of oxygen within the screened area.

Description	Code
Oxygen Delivery Bar & Drape Screen	21365
Clamp - Universal (Socket & Rail) - Handwheel	10250



Push/Pull Bar

This foot-end mounted removable accessory helps to cater for those scenarios where additional manoeuvring requirements arise.

Description	Code
Push/Pull Bar	21390



Theatre Pillow

Tired of the number of pillows heading to the laundry reducing their availability? This superbly comfortable theatre pillow is covered in a water-impermeable material and reduces the turnaround time for pillows between uses. Simply clean & disinfect by hand between uses.

Description	Code
Theatre Pillow	10545




QA3 Accessory Bracket

A new side rail accessory bracket which can be fitted to our QA3™ Stretchers, allowing for other QA3 accessories to be fitted with ease.

Description	Code
QA3 Accessory Bracket	21156

QA3™ OPHTHALMIC STRETCHER

Impressive head-end surgical access, low height, patient comfort and light weight manoeuvrability all combine to produce a perfect transport, treatment and recovery stretcher for ophthalmology.

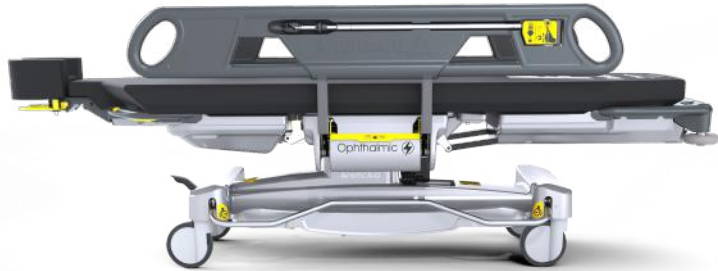

TRANSPORT TREATMENT RECOVERY



Subject to annual servicing




We offer the unique and popular
‘Lifetime Warranty Scheme’

Take out a **AnetiCare** service package for regular and predictable maintenance costs that includes a **10 YEAR WARRANTY**.

Covers engineer’s time and most parts. Contact us to learn more about what an **AnetiCare** service package can offer you.

Available via your Anetic Aid partner – find them at aneticaid.com/distributors.

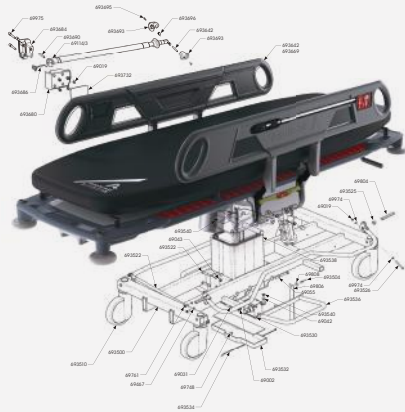

Build options must be specified at the time of ordering



DESIGN & INNOVATION



MADE IN THE UK



GLOBAL DISTRIBUTION



QA4 Mobile Surgery System

The complete multi-purpose solution for transport, treatment and recovery.

PLUS
our other product areas



Scan this code to see our **QA4's** journey through day surgery!



SEE
ALSO



For more information on any of our products or service contracts, call:
sales@aneticaid.com +44 (0)1943 878647 aneticaid.com



DESIGN & INNOVATION



MADE IN THE UK



GLOBAL DISTRIBUTION

Innovative **Medical** Technology - **Practically Applied**



Anetic Aid Ltd, Sapper Jordan Rossi Park, Otley Road, Baildon BD17 7AX, UK