


QA3™

Ophthalmic Stretcher

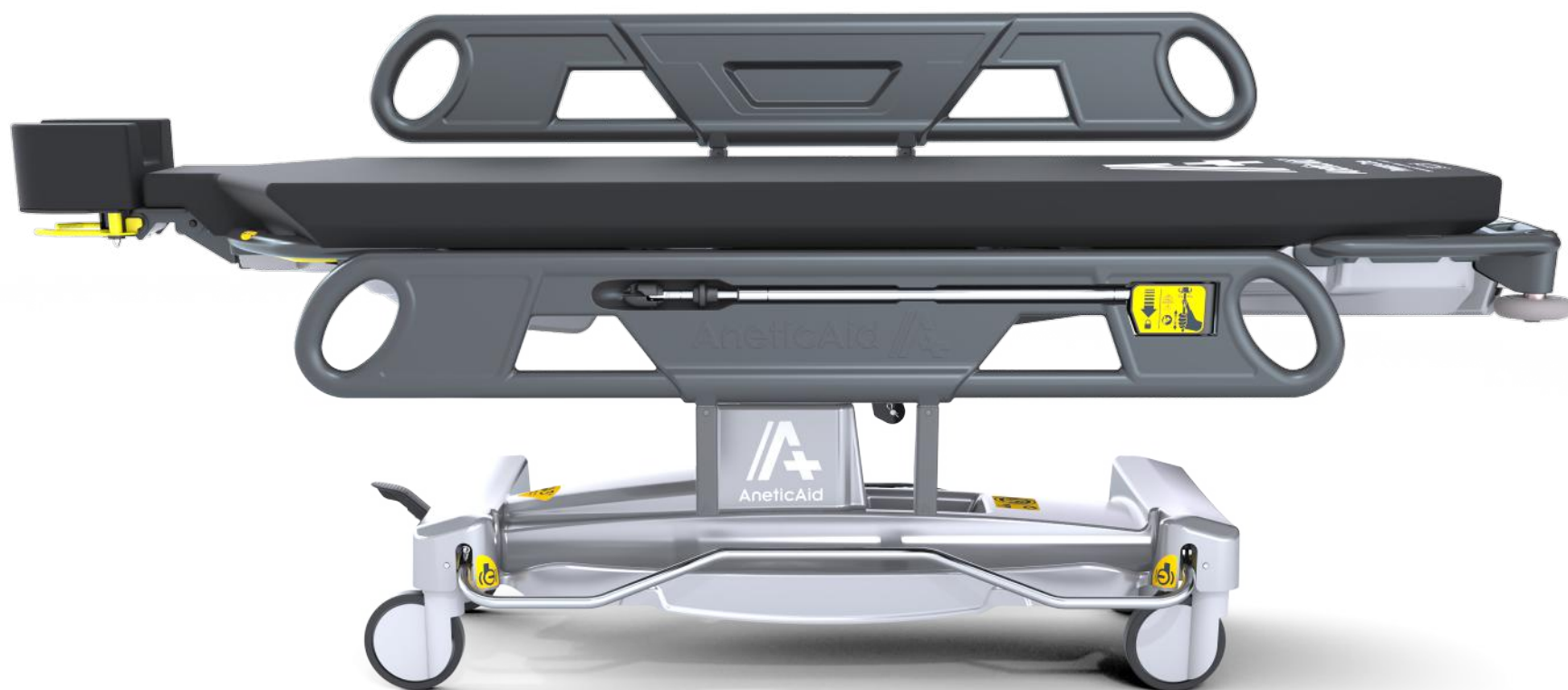

Specification Sheet

AneticAid® 

QA3™

Ophthalmic Stretcher

Dual-articulating headrest

This innovative, intuitive and easy to use head positioner is now fitted as standard, articulating not only on the head but the neck as well. The positioner offers superb surgical positioning, surgical access and patient comfort.

Space saving design

With the back rest and KneeFlex raised, the stretcher has an incredibly small footprint - making it a comfortable and manoeuvrable recovery or transport position and easier to store when not in use.

Fixed transfusion pole

Innovative quick release safety locking mechanism, integrated into the right-hand foot end side rail for easy access and storage.

Six powered functions

Including platform raise and lower, back rest, KneeFlex, Profile, Trendelenburg and Reverse Trendelenburg - as well as CPR return.

For more information on any of our products or service contracts, call:

sales@aneticaid.com **+44 (0)1943 878647** aneticaid.com



DESIGN & INNOVATION



MADE IN THE UK



GLOBAL DISTRIBUTION

Innovative **Medical** Technology - **Practically Applied**

Anetic Aid Ltd, Sapper Jordan Rossi Park, Otley Road, Baildon BD17 7AX, UK

Simultaneous braking low rolling resistance castors

Polyurethane tyres offer low rolling resistance for high manoeuvrability: 150mm diameter central total lock castors fitted as standard, helping to offer a rigid and reliable work platform.

Exceptional patient comfort

Low height of 495mm (excluding mattress) & store away side rails allow for easier mount/dismount, particularly for frail and elderly patients. Also including a versatile positioning headrest, wide patient platform and pressure care mattress for comfort.

Press on / Press off 5th wheel steer

Our unique mechanism ensures maximum manoeuvrability and means the QA3 is easily guided by one support worker.

Pressure care mattress

Manufactured from temperature-sensitive gel-feel foam that moulds to the patient's body shape, relieving pressure by distributing the load as evenly as possible. Outer material is anti-microbial, water resistant, vapour permeable and latex free. Now 5% resistant to chlorine and other disinfectant compounds.



SPECIFICATIONS

Code: 21118

Minimum height	*495mm	Mattress length	2245mm
Maximum height	*875mm	Mattress width	700mm
*excludes mattress		Mattress depth	90mm
Surgical Access	585mm	Maximum patient weight	320kg
Minimum length	1700mm	Backrest	0 to +70°
Maximum length	2275mm	KneeFlex	0 to +30°
Overall width	805mm	Trendelenburg	0 to -12°
Castor diameter	150mm	Reverse Trendelenburg	0 to -10°
Weight	135.5kg	Headrest	0 to -87 / +45°



AneticAid®